


PressureCare™ Mattresses... because YOU care about

- Pressure relief
- Skin protection
- Comfort
- Guaranteed durability

Skil-Care's new **PressureCare™ Mattresses** are value priced and help you provide quality care without compromise. Designed for nursing home and hospital use, **PressureCare™ Mattresses** have these important features...

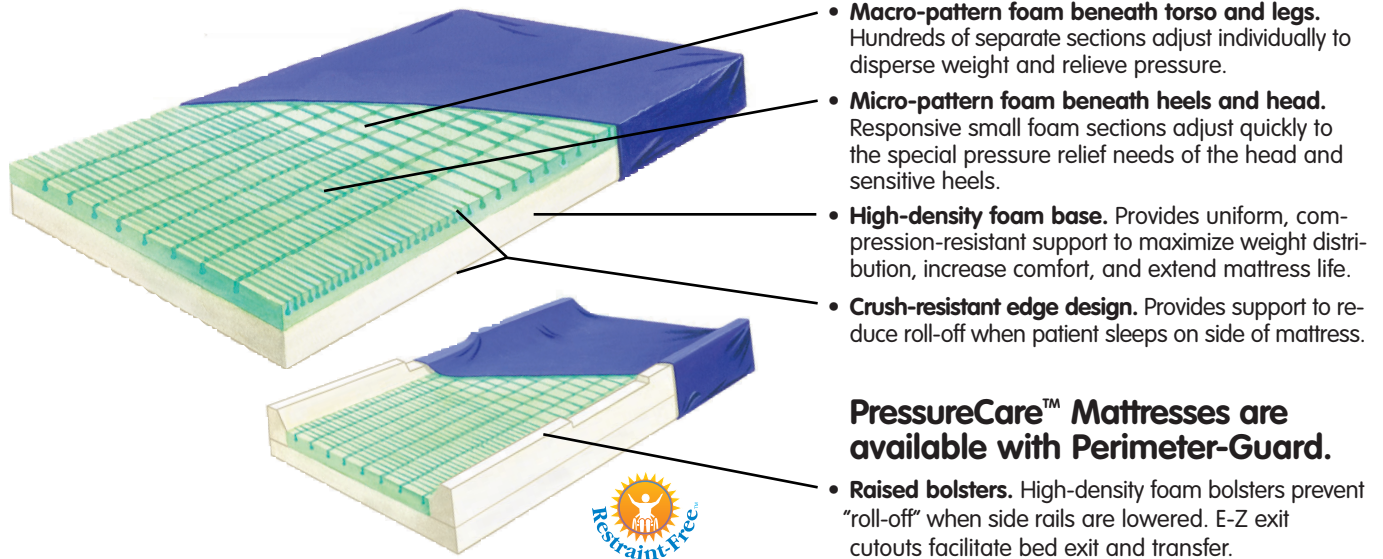

- **Macro-pattern foam beneath torso and legs.** Hundreds of separate sections adjust individually to disperse weight and relieve pressure.
 - **Micro-pattern foam beneath heels and head.** Responsive small foam sections adjust quickly to the special pressure relief needs of the head and sensitive heels.
 - **High-density foam base.** Provides uniform, compression-resistant support to maximize weight distribution, increase comfort, and extend mattress life.
 - **Crush-resistant edge design.** Provides support to reduce roll-off when patient sleeps on side of mattress.
- PressureCare™ Mattresses are available with Perimeter-Guard.**
- **Raised bolsters.** High-density foam bolsters prevent "roll-off" when side rails are lowered. E-Z exit cutouts facilitate bed exit and transfer.



All PressureCare™ Mattresses are available with a choice of two covers:

- **Low-Shear II** A double-stretch, friction-reducing, breathable, fluid-proof fabric that wipes clean for easy care. (Cover guarantee: Two years.)
- **Nylon** Heavy-duty, liquid-proof nylon with polyurethane backing. Wipes clean for easy care. (Cover guarantee: Five years.)

All covers have a slip-resistant bottom for patient safety.

Pressure-Care™ Mattresses comply with fire code 16 CFR 1633.

Seven-year non-prorated guarantee on all foam mattress components and internal fire barrier. Covers are guaranteed separately (see above). Guarantees apply only under conditions of normal use.

Size	Reorder #	Description	Unit	Size	Reorder #	Description	Unit
6" x 36" x 72"	558150	PressureCare	Each	6" x 36" x 72"	558160	Perimeter-Guard	Each
6" x 36" x 76"	558152	PressureCare	Each	6" x 36" x 76"	558162	Perimeter-Guard	Each
6" x 36" x 80"	558154	PressureCare	Each	6" x 36" x 80"	558164	Perimeter-Guard	Each
6" x 36" x 84"	558156	PressureCare	Each	6" x 36" x 84"	558166	Perimeter-Guard	Each

(Mattresses will be shipped with Low-Shear II covers unless otherwise specified. Add the letter "B" to the above model numbers for nylon covers.)

For Information Contact:

Skil-Care™
CORPORATION

29 Wells Ave., Yonkers, NY 10701

Advanced design for nursing and therapy.™

914-963-2040 800-431-2972 Fax 914-963-2567